



dermapose
Refresh

Fat Transfer Technology Dedicated to Volumization & Rejuvenation



www.hcp.dermapose.com


Discover What Is Possible
With Regenerative Fat

Harvest, Purify, Resize and Inject Microfat All Inside Our Closed System Dual-Chamber Syringe



THE QUICKEST OFFICE-BASED PROCEDURE TO OBTAIN QUALITY MICROFAT IN UNDER 60 MINUTES

This sterile and single use syringe allows doctors to reinject ~15cc of microfat. The 50 mL vacuum-lock syringe with a built-in 800 μ m filter designed to gently microsize adipose tissue parcel size allows for easier injection through small 18-21G injectors.

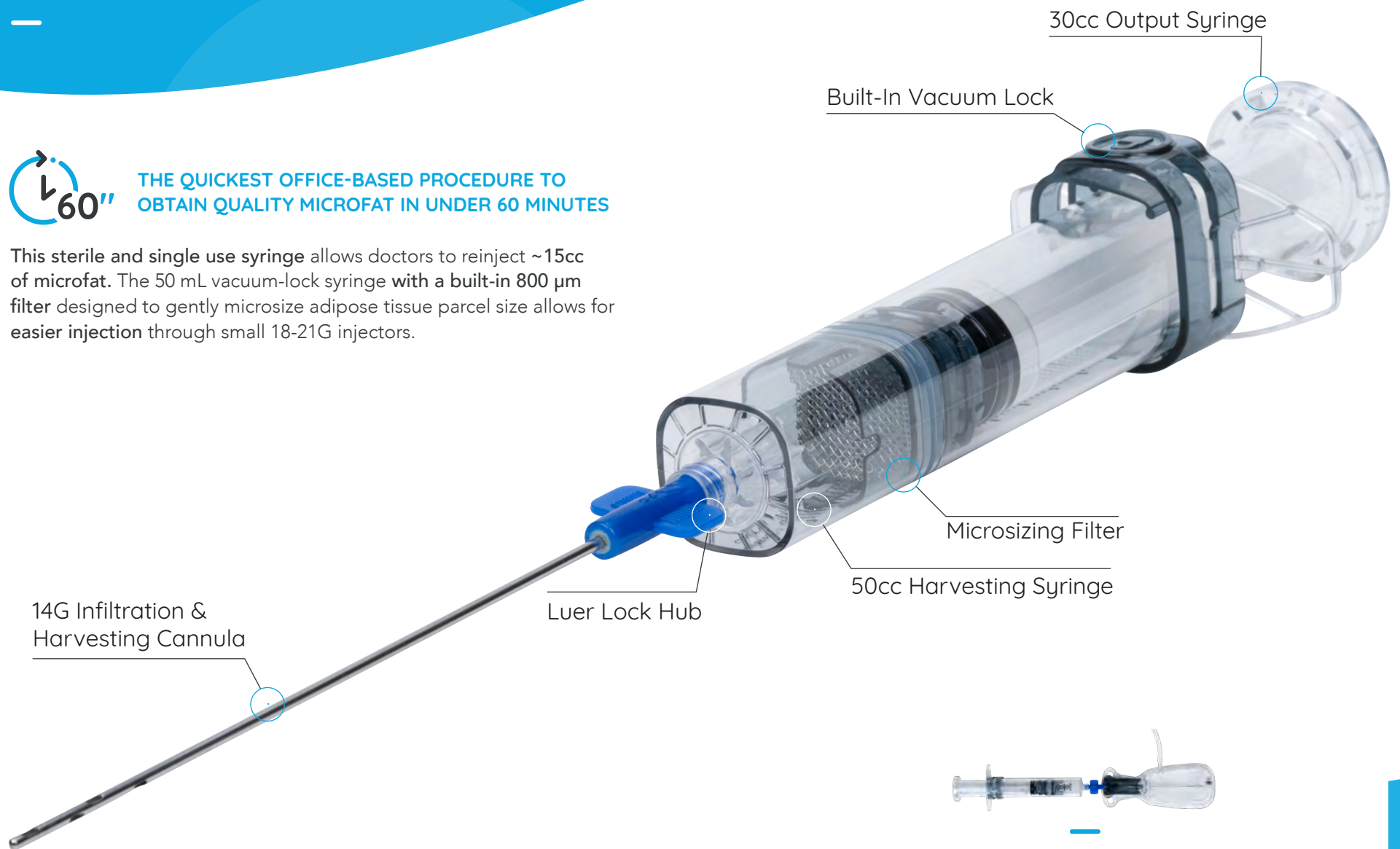

We have also developed a guided harvesting aid that integrates seamlessly with **Dermapose Refresh**.

Discover more on www.hcp.dermapose.com

A closer look at the technology



DUAL MICROSIZING FILTER TECHNOLOGY

Designed to gently step down harvested adipose tissue sizing in a single step for injection through a 21G cannula.

Luer-lock port improves system sterility whilst enabling direct access to purified microfat using the 30cc output syringe.

DERMAPOSE REFRESH STAND

Designed exclusively for Dermalpose Refresh. Provides secure placement of dual chamber syringe during decantation and purification steps. Securely fastens syringe during microsizing.




BIMINI
HEALTHTECH

Bimini Health Tech
420 Stevens Ave, Suite 220
Solana Beach, CA 92075

1-858-348-8050
info@BiminiHealthTech.com
www.BiminiHealthTech.com

The best filtration technology to obtain quality microfat

Designed to microsize tissue in a gentle and reproducible manner, Dermapose preserve the regenerative properties of fat

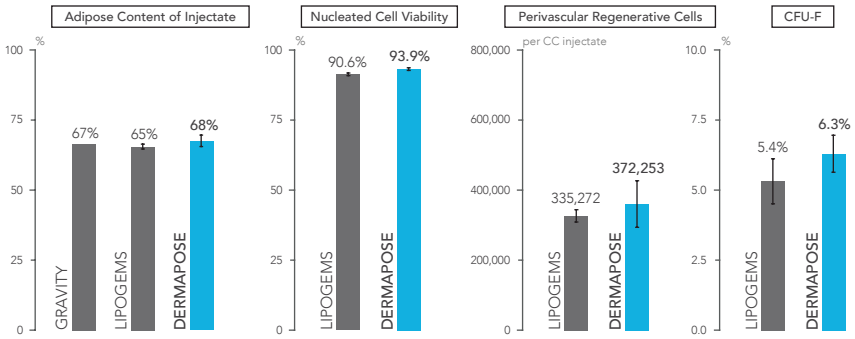

J Magalon, Performed April 20191 at Laboratoire de Culture et Thérapie Cellulaire, Hopital La Conception, Marseille France [Data on file].



“

Being able to **harvest, purify and microsize** 15mL of adipose tissue within a sterile closed-system syringe opens fat grafting up to patients in an office setting.

Dermapose Refresh is the ideal way to obtain microfat in an easy, reproducible manner.

”

E DANIELS, MD

Dermapose Refresh the ideal filler



SIMPLE

Intuitive process enables quick and easy preparation of ~15cc of structural microfat



STERILE

Dual-chamber syringe enables harvesting, purification and microsizing of fat within a closed system, minimizing risk of environmental contaminants



SUPERIOR

Structural microfat sized to 800 microns allowing easy injection through a 21G cannula



dermapose

You need more information?

info@dermapose.com

www.hcp.dermapose.com



BIMINI
HEALTHTECH

Bimini Health Tech
420 Stevens Ave, Suite 220
Solana Beach, CA 92075

1-858-348-8050
info@BiminiHealthTech.com
www.BiminiHealthTech.com

500-039-01-B